

CHROMATE[®] 4300

ELISA 96 well Microwell plate reader



sales@awaretech.com

Features of the Chromate[®]




- ◆ Economical open system & compact design
- ◆ Reads an 96 microwell plate in 12 seconds
- ◆ Multiple calculation modes
- ◆ Dependable performance
- ◆ Custom filter configurations available
- ◆ Intuitive software, designed for ease of use

ChroMate[®] is a powerful yet compact 8-channel microwell plate reader which features: Precision optics, user-programmable assays, standard four filter model (405, 450, 492, and 630 nm) with (optional) six filter VIS model.

ChroMate[®] is an open system and may be used to read most microwell based ELISA kits.

Once installed, ChroMate[®] management software comes preconfigured to interface with your PC and the ChroMate[®] analyzer via a high-speed USB connection.

The ChroMate[®] with your PC provides the rapid reading of 96 well microwell plates or individual microwell strips, customizable reports, calibration curve fit, and unlimited data storage to meet the requirements of the modern laboratory.

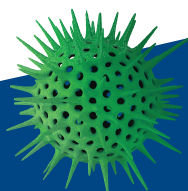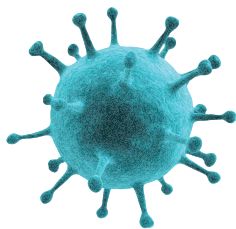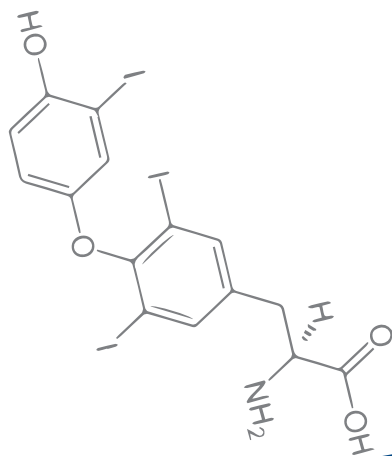

 **AWARENESS[®]**
TECHNOLOGY, INC.

CHROMATE®

4300 Microplate Reader Specifications

Optical measurement range	0.0 to 4.0 absorbance units (A)
Photometric accuracy	± (1% of the reading +/-0.010) from 0 to 1.5 ABS) ± (2% of the reading +/-0.010) from 1.5 to 3.0 ABS)
Stability	Drift of no more than 0.005A in 8 hours
Light source	Halogen lamp source with lamp saver feature
Wavelengths	Standard: 405, 450, 492, and 630 nm Optional 6-filter: 405, 450, 492, 545, 600, and 630 nm (must be ordered at time of manufacturing) Filter options from 405 nm to 750 nm
Filter Type:	IAD hardcoat interference, 10 nm half bandpass
Vessel:	Standard 96-well microwell plates, or strip trays
Power Requirements:	100V or 240V AC, 70W, 50-60Hz (auto-sensing)
PC connection:	USB Port
Speed:	Reads absorbances of 96 wells in about 12 seconds
Plate Transport:	Stepper motor
PC interface with Chromate	Manager software - included. (PC not included)
Calculation Modes:	Single point calibration, point-to-point curve fit, polynomial regression, linear and sigmoidal regression (log & linear), cutoff absorbance, multi-point % absorbance
Additional Features:	User-programmable open system, selectable plate formatting, alpha numeric test naming, automatic interpretation options, duplicate well options, curve plotting and editing, flags and error messages
Certifications:	NRTL Listed, CE Mark Awareness Technology Quality System is certified under ISO13485:2016
Enclosure:	All metal enclosure
Dimensions:	Approx. 7.5x8.0x16.5 (19x20.3x42cm) Weights 11.5 lbs. (5.22 kg)

ISO 13485:2016



MADE IN U.S.A.

Dri-Dye® Check Strips

For performance verification of Stat Fax®, Chromate®, and other microstrip and microplate readers



Awareness Technology, Inc. USA & International
tel: +1 772 283 6540 fax: +1 772 283 8020

Awareness Technology, Inc. Europe
tel: +43 (0) 2236 892465 fax: +43 (0) 2236 892464

Awareness Technology, Inc. Asia & Africa
tel: +971 (06) 557 80 58 fax: +971 (06) 557 80 59

